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MAX	Before use turn switch to max position. (C)
MIN	After use turn switch to min position (C)
Lampetest	Control the alarm lamp by pushing the lampetest button for 3-4 sec. The lamps should be lit, if not notify. (C)
Red flashing light	Indicate to low airflow. Indicate to low airflow or to high front position. – CLOSE. Notify. (A)
Correct sash opening	Correct opening of sash is indicated with stickers. (B)

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
Drift	The fume hood is always on.
Grønt lys	Correct air-flow
Max	Maximum airflow regardless of position of front.
Red blink + lyd	Indicate to low airflow. Indicate to low airflow or to high front position. – CLOSE
Red lys + lyd	Red lamp glows constantly, sound for 30 secs; electricity problem – CLOSE
Alarm	Turns alarm on or off.
Light	Turns light on or off.
<i>All LABFLEX Fume hoods are controlled by LAB-VENT control panel</i>	

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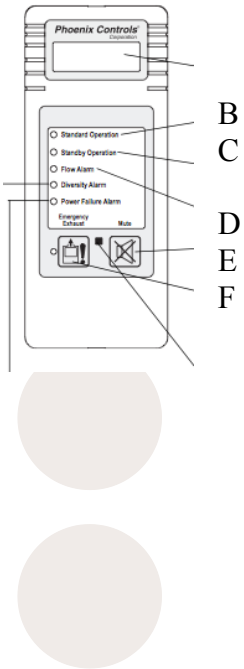Standard Operation LED	When this LED is on, it indicates the hood is operating at the standard face velocity or flow. (B)
Standby Operation LED	When this LED is on, it indicates the hood is operating at a lower face velocity or flow (active only with a ZPS). (C)
Flow Alarm/Drive Failure LED	When this LED is on, it indicates unsafe airflow condition. (D)
Mute button	Push the button to silence an alarm. The Flow Alarm LED will remain on. Mute mode is reset when the alarm conditions clear. (E)
Emergency Exhaust button with LED—	Push this button to activate the emergency exhaust mode. The alarm will sound and the LED at the left of the Emergency Exhaust button will blink. In this mode, exhaust air is at its maximum flow. Push the button again to turn off the emergency exhaust/test mode. Note: On constant volume hoods, this button is used only to test the alarm circuit, not modulate the exhaust valve. (F)

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Fume hood controlled by "Phoenix controls"

Face Velocity or Flow Display Indicates the face velocity or flow of the hood. (Available only on the FHM631.) (A)

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
Flow Alarm/Drive Failure LED When this LED is on, it indicates unsafe airflow condition. (D)

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
Nødventilasjon Push this button to activate the emergency exhaust mode. The alarm will sound and the LED at the left of the Emergency Exhaust button will blink. In this mode, exhaust air is at its maximum flow. Push the button again to turn off the emergency exhaust/test mode.
Note: On constant volume hoods, this button is used only to test the alarm circuit, not modulate the exhaust valve. (F)

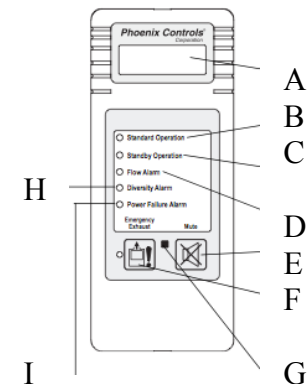Light sensor	Detects the light level in the room. Triggers the energy waste alert, if configured. (Available only on the FHM631.) (G)
Diversity Alarm LED	Alerts lab users to reduce the total flow by closing their sashes. This visual alarm is triggered when the flow demand exceeds the flow limit and a diversity alarm is generated by the Celeris system or the BMS. (H)
Power Failure Alarm LED	Activates when there is a loss of power. Used in conjunction with the Power Failure Alarm. With this option, a separate PLA module is provided with the monitor. (I)

Face Velocity or Flow Display	Indicates the face velocity or flow of the hood. (Available only on the FHM631.) (A)
Standard Operation LED	When this LED is on, it indicates the hood is operating at the standard face velocity or flow. (B)
Standby Operation LED	When this LED is on, it indicates the hood is operating at a lower face velocity or flow (active only with a ZPS). (C)

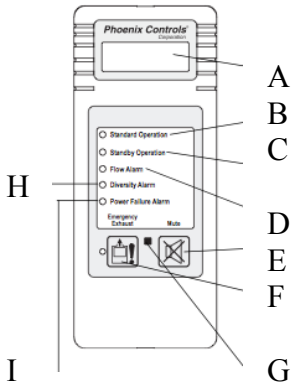Flow Alarm/Drive Failure LED	When this LED is on, it indicates unsafe airflow condition. (D)
Mute button	Push the button to silence an alarm. The Flow Alarm LED will remain on. Mute mode is reset when the alarm conditions clear. (E)
Emergency Exhaust button with LED—	Push this button to activate the emergency exhaust mode. The alarm will sound and the LED at the left of the Emergency Exhaust button will blink. In this mode, exhaust air is at its maximum flow. Push the button again to turn off the emergency exhaust/test mode. Note: <i>On constant volume hoods, this button is used only to test the alarm circuit, not modulate the exhaust valve. (F)</i>
Light sensor—	Detects the light level in the room. Triggers the energy waste alert, if configured. (Available only on the FHM631.) (G)
Diversity Alarm LED—	Alerts lab users to reduce the total flow by closing their sashes. This visual alarm is triggered when the flow demand exceeds the flow limit and a diversity alarm is generated by the Celeris system or the BMS. (H)
Power Failure	Activates when there is a loss of power. Used in conjunction with the Power Failure Alarm. With this

Alarm LED option, a separate PLA module is provided with the monitor. (I)

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General rules:

1. Maximum working opening is indicated
2. Always use as low opening as possible.
3. Avoid large set up. This disrupts the ventilation flow
4. Place work set up as close as possible to the back - never along side of the cabinet
5. Keep cabinet clean and tidy
6. No not use for storage
7. Use suitable clothing
8. Is failures detected, stop working immediately. "Out of order" sign is placed and responsible person is contacted.